

SENTIFIT® pierceTube



For the quantitative detection of human hemoglobin in stool samples

Instructions


The SENTIFIT® pierce Tube is extremely easy and hygienic to use. The following instruction shows how to collect a stool sample correctly in 6 steps.

You can find a short video on stool collection at:

[Darmkrebsscreening](#)



Step 1

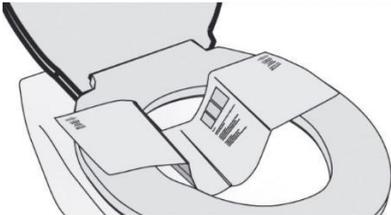


Attach the stool catcher to the toilet seat. Defecate over the paper of the stool catcher. Avoid contact with water or urine.

Step 2



Unscrew the green cap containing the applicator stick. Do not hold the tube upside down.

Step 3

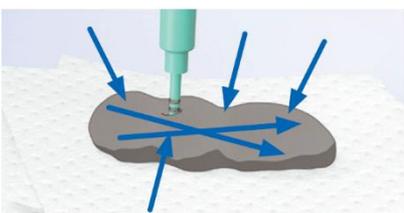

Insert the sample stick into the stool at four different sites, being careful not to go too deep. Then swipe the stick across the surface in a crisscross pattern.

SENTIFIT®pierceTube – Soeinfach geht's

Step 4



Make sure that all the notches are filled with stool.
Do not use too much stool.

Step 5



Screw the green cap back onto the tube.

Step 6



Tilt the tube gently upside down several times.
Do not open it after this point. Return
the tube to your doctor no later than the day after
the sample is taken.

General Information

- Keep the container out of the reach of children.
- The liquid in the container must not be swallowed and must not come into contact with the eyes or mucous membranes.
- The container is designed to prevent the liquid from leaking.
- In case of accidental contact, seek medical advice immediately.
- Collect the stool sample immediately after defecation.
- Do not insert the container or any part of it into the anus.
- Stool samples must not be collected during menstruation, if there is bleeding from the anus (e.g. due to hemorrhoids or other causes), or if the stool is contaminated with urine.
- Use the container only once and avoid direct sunlight.
- Do not open the transparent cap.